For Research Use Only

CTSO Polyclonal antibody Catalog Number:14296-1-AP


Antibodies | ELISA kits | Proteins WWW.ptglab.com

Basic Information	Catalog Number: 14296-1-AP	GenBank Accession Number: BC049206		Purification Method: Antigen affinity purification	
	Size: 150ul, Concentration: 450 ug/ml by Nanodrop and 260 ug/ml by Bradford method using BSA as the standard; Source: Rabbit	GenelD (NCBI): 1519		Recommended Dilutions: IP 0.5-4.0 ug for 1.0-3.0 mg of total	
		UNIPROT ID: P43234		protein lysate IHC 1:50-1:500	
		Full Name: cathepsin O			
	Isotype: IgG	Calculated MW: 36 kDa			
	Immunogen Catalog Number: AG5651	Observed MW: 36 kDa			
Applications	Tested Applications: IP, IHC, ELISA		Positive Controls: IP : HEK-293 cells,		
	Species Specificity: human, mouse, rat		IHC : human ovary tumor tissue,		
	Note-IHC: suggested antigen retrieval with TE buffer pH 9.0; (*) Alternatively, antigen retrieval may be performed with citrate buffer pH 6.0				
Storage	Storage: Store at -20°C. Stable for one year after shipment. Storage Buffer: PBS with 0.02% sodium azide and 50% glycerol pH 7.3.				
*** 20ul sizes contain 0.1% BSA	Aliquoting is unnecessary for -20 $^\circ$ C s	torage			


For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA)E: proteintech@ptglab.comW: ptglab.comW: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


Selected Validation Data

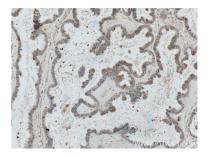

Immunohistochemical analysis of paraffinembedded human ovary tumor tissue slide using 14296-1-AP (CTSO antibody) at dilution of 1:200 (under 10x lens. Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0).

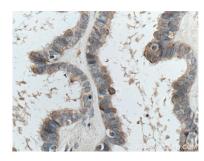

Immunohistochemical analysis of paraffinembedded human ovary tumor tissue slide using 14296-1-AP (CTSO antibody) at dilution of 1:200 (under 40x lens. Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0).

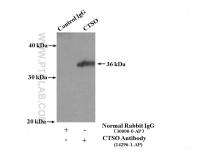

IP result of anti-CTSO (IP:14296-1-AP, 4ug; Detection:14296-1-AP 1:300) with HEK-293 cells lysate 3600ug.