For Research Use Only

CoraLite® Plus 647-conjugated DKK3 Polyclonal antibody



Purification Method:

Catalog Number: CL647-10365

Featured Product

Basic Information

Catalog Number: GenBank Accession Number:

CL647-10365 BC007660 Antigen affinity purification GeneID (NCBI): Recommended Dilutions: 100ul , Concentration: 1000 $\mu g/ml$ by 27122 IF/ICC 1:50-1:500

Nanodrop: **UNIPROT ID:** Excitation/Emission maxima

Q9UBP4 wavelengths: 654 nm / 674 nm Rabbit Full Name:

Isotype: dickkopf homolog 3 (Xenopus laevis)

IgG Calculated MW:

Immunogen Catalog Number: 38 kDa AG0417

Observed MW:

55 kDa

Applications

Tested Applications: IF/ICC, FC (Intra) Species Specificity:

human, mouse

Positive Controls:

IF/ICC: SH-SY5Y cells,

Background Information

DKK3, also named as REIC, belongs to the dickkopf family. It antagonizes canonical Wnt signaling by inhibiting LRP5/6 interaction with Wnt and by forming a ternary complex with the transmembrane protein KREMEN that promotes internalization of LRP5/6. DKKs play an important role in vertebrate development, where they locally inhibit Wnt regulated processes such as antero-posterior axial patterning, limb development, somitogenesis and eye formation. In the adult, Dkks are implicated in bone formation and bone disease, cancer and Alzheimer's disea. The adenovirus carrying REIC/Dkk-3 (Ad-REIC) exhibits a potent tumor-specific cell-killing function for various human cancers. It has also become evident that some human cancers are resistant to Ad-REIC-induced apoptosis. The DKK3 is a Glycosylation protein with MW about 55 kDa or 38-43 kDa.

Storage


Store at -20°C. Avoid exposure to light. Stable for one year after shipment. Storage Buffer

PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.


Aliquoting is unnecessary for -20°C storage

in USA), or 1(312) 455-8498 (outside USA)

Selected Validation Data

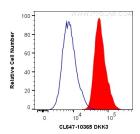

1X10^6 HepG2 cells were intracellularly stained with 0.2 ug CoraLite® Plus 647 Anti-Human DKK3 (CL647-10365) (red), or 0.2 ug isotype control (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).

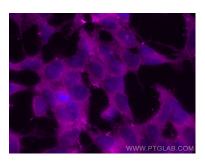

Immunofluorescent analysis of (-20°C Ethanol) fixed SH-SY5Y cells using CoraLite® Plus 647 DKK3 antibody (CL647-10365) at dilution of 1:200.