

## Human TNFRSF9/CD137 Sandwich ELISA Kit Datasheet

Please read it entirely before use

**Catalogue Number:** KE00283

**Size:** 96T

**Sensitivity:** 4.4 pg/mL

**Range:** 15.6-1000 pg/mL

**Usage:** For the quantitative detection of human TNFRSF9/CD137 concentrations in serum, plasma and cell culture supernatant.


This product is for research use only and not for use in human or animal therapeutic or diagnostic.

Table of content	page
1. Background	3
2. Principle	3
3. Required Materials	3
4. Kit Components and Storage	4
5. Safety Notes	4
6. Sample Collection and Storage	4
7. Regent Preparation	5
8. Assay Procedure Summary	6
9. Validation Data	7
9.1 Standard curve	7
9.2 Precision	7
9.3 Recovery	8
9.4 Sample values	8
9.5 Sensitivity	8
9.6 Linearity	9
9.7 Specificity	9
10. References	10

## 1. Background

CD137, also known as TNFRSF9 or 4-1BB, is an inducible T cell surface receptor which belongs to the tumor necrosis factor receptor superfamily. CD137 is a transmembrane protein expressed on the surface of activated T-cells. In addition, activation-dependent expression of CD137 has also been found in B lymphocytes, monocytes, and diverse nonlymphoid cell types. CD137 provides a co-stimulatory signal that enhances the survival, and differentiation of cells, and has a crucial role in the development of CD8 cytotoxic T cells and anti-tumor immunity. Soluble forms of CD137 (sCD137) are generated by differential splicing and are released by activated lymphocytes. Elevated serum sCD137 level has been associated with some diseases, including rheumatoid arthritis, multiple sclerosis, chronic lymphocytic leukemia, and acute coronary syndrome.

## 2. Principle

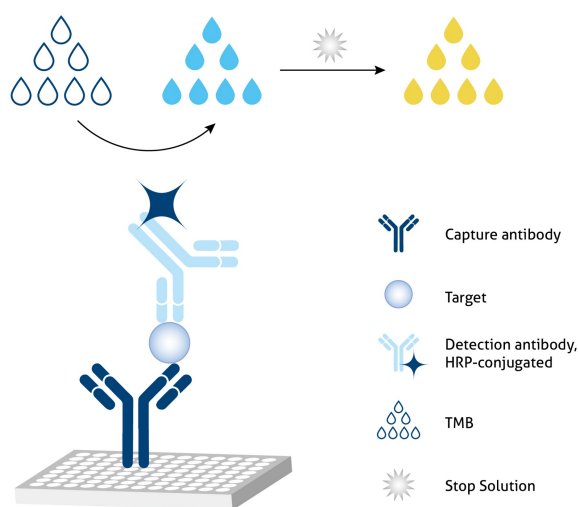

### Sandwich ELISA structure (Detection antibody labeled with HRP)

A capture antibody is pre-coated onto the bottom of wells which binds to analyte of interest. A detection antibody labeled with HRP also binds to the analyte. TMB acts as the HRP substrate and the solution color will change from colorless to blue. A stop solution containing sulfuric acid turns solution yellow. The color intensity is proportional to the quantity of bound protein which is measurable at 450 nm with the correction wavelength set at 630 nm.

## 3. Required Materials

- 3.1 A microplate reader capable of measuring absorbance at 450 nm with the correction wavelength set at 630 nm.
- 3.2 Calibrated, adjustable precision pipettes and disposable plastic tips. A manifold multi-channel pipette is recommended for large assays.
- 3.3 Plate washer: automated or manual.
- 3.4 Absorbent paper towels.
- 3.5 Glass or plastic tubes to prepare standard and sample dilutions.
- 3.6 Beakers and graduated cylinders.
- 3.7 Log-log or semi-log graph paper or computer and software for ELISA data analysis. A four-parameter logistic (4-PL) curve-fit is recommended.

## 4. Kit Components and Storage

Microplate - antibody coated 96-well microplate (8 well × 12 strips)	1 plate	<b>Unopened Kit:</b>  Store at 2-8°C for 6 months or -20°C for 12 months.  <b>Opened Kit:</b>  All reagents stored at 2-8°C for 7 days.  <b>Please use a new standard for each assay.</b>
Protein standard - 2000 pg/bottle; lyophilized	2 bottles	
Detection antibody, HRP-conjugated (100×) - 120 µL/vial*	1 vial	
Sample Diluent PT 1 - 30 mL/bottle	1 bottle	
Detection Diluent - 30 mL/bottle	1 bottle	
Wash Buffer Concentrate (20×) - 30 mL/bottle	1 bottle	
Tetramethylbenzidine Substrate (TMB) - 12 mL/bottle	1 bottle	
Stop Solution - 12 mL/bottle	1 bottle	
Plate Cover Seals	4 pieces	

## 5. Safety Notes

- 5.1 Avoid any skin and eye contact with Stop Solution and TMB. In case of contact, wash thoroughly with water.
- 5.2 Do not use the kit after the expiration date.
- 5.3 Do not mix or substitute reagents or materials from other kit lots or other sources.
- 5.4 Be sure to wear protective equipment such as gloves, masks and goggles during the experiment.
- 5.5 When using an automated plate washer, adding a 30 second soak period following the addition of Wash Buffer to improve assay precision

## 6. Sample Collection and Storage

- 6.1 Serum: Allow blood samples to clot for 30 minutes, followed by centrifugation for 15 minutes at 1000xg. Clear serum can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.
- 6.2 Plasma: Use EDTA, heparin, or citrate as an anticoagulant for plasma collection. Centrifuge for 15 minutes at 1000xg within 30 minutes of collection. The plasma can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.
- 6.3 Cell Culture Supernatant: Remove particulates by centrifugation for 5 minutes at 500xg and assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

## 7. Regent Preparation

**7.1 Wash Buffer (1X):** If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Add 30 mL of Wash Buffer Concentrate(20X) to 570 mL deionized or distilled water to prepare 1X Wash Buffer.


**7.2 Detection Antibody, HRP-conjugated(1X):** Dilute 100X Detection Antibody, HRP-conjugated 1:100 using Detection Diluent prior to assay. Suggested 1:100 dilution: 10  $\mu$ L 100X Detection Antibody, HRP-conjugated + 990  $\mu$ L Detection Diluent (Centrifuge the 100 X Detection Antibody solution, HRP-conjugated for a few seconds prior to use)

**7.3 Sample Dilution:** Different samples should be diluted with corresponding Sample Diluent, samples may require further dilution if the readout values are higher than the highest standard OD reading. Variations in sample collection, processing and storage may affect the results of the measurement.

Recommended Dilution for different sample types: 1:2 or 1:4 is recommended for human serum and plasma; 1:2 or 1:4 is recommended for cell culture supernatant.

### 7.4 Standard Serial Dilution:

Add 2 mL Sample Diluent PT 1 in protein standard.

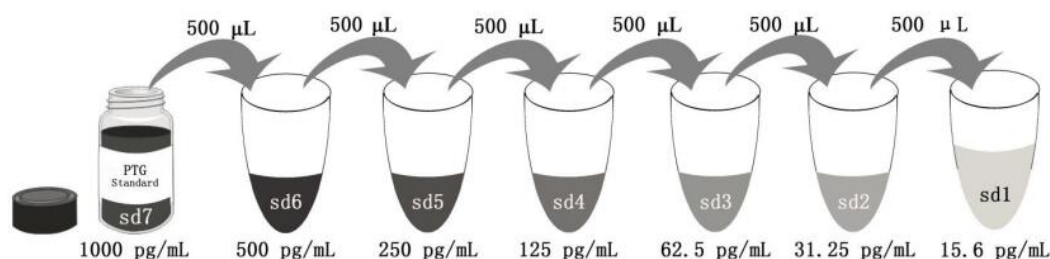

Add # $\mu$ L of Standard diluted in the previous step	—	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L
# $\mu$ L of Sample Diluent PT 1	2000 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L
	"sd7"	"sd6"	"sd5"	"sd4"	"sd3"	"sd2"	"sd1"

## 8. Assay Procedure Summary

Bring all reagents to room temperature before use (Detection antibody, HRP-conjugated can be used immediately). To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

8.1 Take out the required number of microplate strips and return excess strips to the foil pouch containing the drying reagent pack and reseal; store at 4°C immediately. Microplate strips should be used in one week.

8.2 Preset the layout of the microplate, including control group, standard group and sample group, add 100 µL of each standard and sample to the appropriate wells. (Make sure sample addition is uninterrupted and completed within 5 to 10 minutes, It is recommended to assay all standards, controls, and samples in duplicate).

8.3 Seal plate with cover seal, pressing it firmly onto top of microwells. Incubate the plate for 120 minutes at 37°C.

### 8.4 Wash

1) Gently remove the cover seal. Discard the liquid from wells by aspirating or decanting. Remove any residual solution by tapping the plate a few times on fresh paper towels.

2) Wash 4 times with 1X Wash Buffer, using at least 350-400 µL per well. Following the last wash, firmly tap plates on fresh towels 10 times to remove residual Wash Buffer. Avoid getting any towel fibers in the wells or wells drying out completely.

8.5 Add 100 µL of 1X Detection antibody, HRP-conjugated solution (refer to Reagent Preparation 7.2) to each well. Seal plate with cover seal and incubate for 40 minutes at 37°C.

8.6 Repeat wash step in 8.4.

8.7 Signal development: Add 100 µL of TMB substrate solution to each well, protected from light. Incubate for 15 to 20 minutes. Substrate Solution should remain colorless until added to the plate.

8.8 Quenching color development: Add 100 µL of Stop Solution to each well in the same order as addition of the TMB substrate. Mix by tapping the side of the plate gently. NB: Avoid skin and eye contact with the Stop solution.

8.9 Read results: Immediately after adding Stop solution read the absorbance on a microplate reader at a wavelength of 450 nm. If possible, perform a double wavelength readout (450 nm and 630 nm).


8.10 Data analysis: Calculate the average of the duplicate readings (OD value) for each standard and sample, and subtract the average of the zero standard absorbance. Construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis, use four-parameter logistic curve-fit (4-PL) analysis to do this. If the samples have been diluted, the OD readout from the standard curve must be multiplied by the dilution factor used.

Step	Reagent	Volume	Incubation	Wash	Notes
1	Standard and Samples	100 µL	120 min	4 times	Cover Wells incubate at 37°C
2	Diluent Detection antibody, HRP-conjugated Solution	100 µL	40 min	4 times	Cover Wells incubate at 37°C
3	TMB Substrate	100 µL	15-20 min	Do not wash	Incubate in the dark at 37°C
4	Stop Solution	100 µL	0 min	Do not wash	-
5	Read plate at 450 nm and 630 nm immediately after adding Stop solution. DO NOT exceed 5 minutes.				

## 9. Validation Data

### 9.1 Standard curve

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

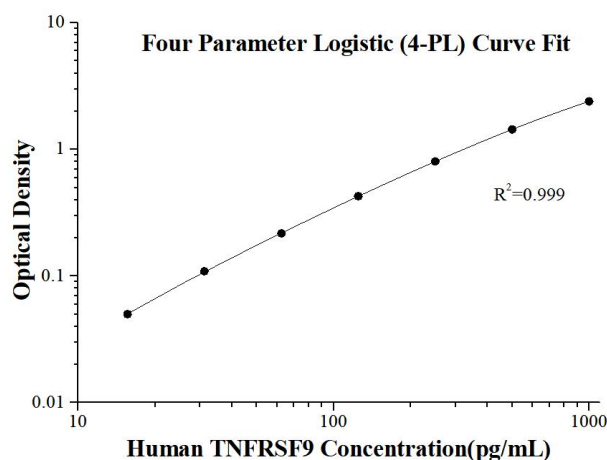

(pg/mL)	O.D	Average	Corrected
0	0.087 0.087	0.087	-
15.6	0.136 0.138	0.137	0.050
31.2	0.201 0.191	0.196	0.109
62.5	0.316 0.291	0.304	0.217
125	0.517 0.511	0.514	0.427
250	0.917 0.865	0.891	0.804
500	1.556 1.496	1.526	1.439
1000	2.545 2.422	2.484	2.397

### 9.2 Precision

**Intra-assay Precision** (Precision within an assay) Three samples of known concentration were tested 20 times on one plate to assess intra-assay precision.

**Inter-assay Precision** (Precision between assays) Three samples of known concentration were tested in 24 separate assays to assess inter-assay precision.

Intra-assay Precision				
Sample	n	Mean (pg/mL)	SD	CV%
1	20	422.9	20.4	4.8
2	20	108.0	3.6	3.3
3	20	28.2	2.1	7.6

Inter-assay Precision				
Sample	n	Mean (pg/mL)	SD	CV%
1	24	494.1	36.3	7.3
2	24	116.8	8.8	7.5
3	24	31.9	3.1	9.7

## 9.3 Recovery

The recovery of human TNFRSF9/CD137 spiked to three different levels throughout the range of the assay in various matrices was evaluated.

Sample Type		Average% of Expected	Range (%)
Human serum	1:2	92	80-106
	1:4	99	91-109
Cell culture supernatant	1:2	103	84-121
	1:4	104	94-120

## 9.4 Sample values

**Serum** - Human serum samples were evaluated for the presence of human TNFRSF9/CD137 in this assay.

Sample Type	Mean (pg/mL)	Range (pg/mL)
Human serum (n=16)	182.9	77.0-548.0

**Cell culture supernatant** - Human peripheral blood mononuclear cells (PBMC) ( $1 \times 10^6$  cells/mL) were cultured in RPMI-1640 supplemented with 10% fetal bovine serum, 100 U/mL penicillin and 100 µg/mL streptomycin sulfate, the cell culture supernatant was stimulated with 10ug/mL of PHA for 3 day and 5 day. An aliquot of the culture supernatant was removed and assayed for levels of human TNFRSF9/CD137.

Condition	Day 3 (pg/mL)	Day 5 (pg/mL)
Unstimulated	32.0	23.2
Stimulated	450.8	527.2

## 9.5 Sensitivity

The minimum detectable dose of human TNFRSF9/CD137 is 4.4 pg/mL. This was determined by adding two standard deviations to the concentration corresponding to the mean O.D. of 20 zero standard replicates.



## 9.6 Linearity

To assess the linearity of the assay, samples were diluted with the appropriate **Sample Diluent** to produce samples with values within the dynamic range of the assay.

		Human serum	Cell culture supernatant
1:2	Average% of Expected	100	100
	Range (%)	-	-
1:4	Average% of Expected	105	107
	Range (%)	96-118	106-109
1:8	Average% of Expected	106	109
	Range (%)	94-117	102-118
1:16	Average% of Expected	107	118
	Range (%)	90-116	114-121

## 9.7 Specificity

This assay recognizes natural and recombinant human TNFRSF9/CD137.

The following factors prepared at 50 ng/mL were assayed and exhibited no cross-reactivity or interference.

Recombinant human:	Recombinant mouse
4-1BB	TNF RI
	TNF RII

A sample containing 50 ng/mL of the recombinant human TNF RI reads as 25.7 pg/mL (0.05% cross-reactivity).

A sample containing 50 ng/mL of the recombinant human TNF RII reads as 599.2 pg/mL (1.2% cross-reactivity)

## 10. References

1. Vinay DS, et al. Role of 4-1BB in immune responses. *Semin Immunol.* 10(6):481-9 (1998).
2. Schwarz H, et al. ILA, the human 4-1BB homologue, is inducible in lymphoid and other cell lineages. *Blood.* 85(4):1043-52 (1995).
3. Nakaima Y, et al. CD137 is induced by the CD40 signal on chronic lymphocytic leukemia B cells and transduces the survival signal via NF- $\kappa$ B activation. *PLoS One.* 16;8(5):e64425 (2013).
4. Michel J, et al. A soluble form of CD137 (ILA/4-1BB), a member of the TNF receptor family, is released by activated lymphocytes and is detectable in sera of patients with rheumatoid arthritis. *Eur J Immunol.* 28(1):290-5 (1998).
5. Sharief MK, et al. Heightened intrathecal release of soluble CD137 in patients with multiple sclerosis. *Eur J Neurol.* 9(1):49-54 (2002).
6. Furtner M, et al. Levels of soluble CD137 are enhanced in sera of leukemia and lymphoma patients and are strongly associated with chronic lymphocytic leukemia. *Leukemia.* 19(5):883-5 (2005).
7. Yan J, et al. Clinical implications of elevated serum soluble CD137 levels in patients with acute coronary syndrome. *Clinics (Sao Paulo).* 68(2):193-8 (2013).