SARS-CoV-2 Spike Rekombinanter Antikörper

Katalog-Nr.:91339-PTG




www.ptglab.com

Allgemeine Informationen	Katalog-Nr.: 91339-PTG	GenBank-Zugangsnummer: NC_045512	Reinigungsmethode: Protein-A-Chromatographie
	Größe: 100ug	GenelD (NCBI): 43740568	CloneNo.: AM038105
	Wirt: Mensch Isotyp: IgG1	Vollständiger Name: SARS-CoV-2 Spike Protein Berechneté Masse: 141 kDa	Empfohlene Verdünnungen: Sample dependent. To be determined by the end user.
Getestete Reaktivität: Virus			
Hintergrundinformationen	COVID-19, which is short for coronavirus disease 2019, is the official name of the respiratory disease caused by infection with the novel coronavirus SARS-CoV-2. The virus that causes COVID-19 was named SARS-CoV-2 because it is a coronavirus genetically similar to, yet distinct from, the virus that caused the severe acute respiratory syndrome (SARS) outbreak in 2003. Studying the details of how this virus replicates and causes the disease will allow scientists and physicians to more rapidly develop fast and accurate methods of detection as well as to deploy therapeutic and vaccine strategies. This antibody was derived from COVID-19 patients who have cleared the virus. Patient serum IgG was sequenced and expressed as full-length IgG1 with human immunoglobulin heavy and light chains in mammalian 293 cells.		
Lagerung	Lagerungsbedingungen: Avoid repeated freeze/thaw cycles by aliquoting items into single-use fractions for storage at -20°C for up to 2 years. Keep all reagents on ice when not in storage. This product is guaranteed for 12 months from date of receipt. Lagerungspuffer: 140 mM Hepes, pH 7,5, 70 mM NaCl, 32 mM NaOAc, 0,035% Natriumazid und 30% Glyzerin.		
*** 20ul-Größen enthalten 0.1% BSA	Aliquotieren ist nicht notwen	dig bei -20 [°] C lagerung	

For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA)E: proteintech@ptglab.comW: ptglab.comW: ptglab.com

This product is made by Active Motif® and available through Proteintech Group, Inc.

Ausgewählte Validierungsdaten

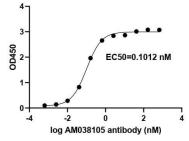

SARS-CoV-2 Spike Antibody (clone AM038105) tested by ELISA.SARS-CoV-2 Spike RBD protein was coated onto microtiter plates at 0.5 µg/mL and then incubated with a dilution series of SARS-CoV-2 Spike Antibody (clone AM009105). Bound antibodies were detected with anti-human IgG conjugated to horseradish peroxidase (HRP) followed by incubation with HRP Substrate and then measuring the resulting absorbance at 450 nm. Data provided by Active Motif®..