

Nur für Forschungszwecke

VDAC1/Porin Polyklonaler Antikörper



Katalog-Nr.:CL647-55259

Vorgestelltes Produkt

Allgemeine Informationen

Katalog-Nr.:

CL647-55259

Größe:

100ul , Konzentration: 1000 µg/ml von7416

Nanodrop;

Wirt:

Kaninchen

Isotyp:

IgG

GenBank-Zugangsnummer:

NM_003374

GeneID (NCBI):

1000

Vollständiger Name:

voltage-dependent anion channel 1

Berechnete Masse:

31 kDa

Beobachtete Masse:

31 kDa

Reinigungsmethode:

Antigen-Affinitätsreinigung

Anregungs-/Emissionsmaxima-

Wellenlängen:

654 nm / 674 nm

Anwendungen

Geprüfte Anwendungen:

FC (Intra)

Getestete Reaktivität:

Human, Maus, Ratte

Hintergrundinformationen

VDAC1, also named as VDAC, porin 31HM, porin 31HL and plasmalemmal porin, belongs to the eukaryotic mitochondrial porin family. It adopts an open conformation at low or zero membrane potential and a closed conformation at potentials above 30-40 mV, to form a channel through the mitochondrial outer membrane and also the plasma membrane. Unlike other membrane transport proteins, porins are large enough to allow passive diffusion. Studies have shown that VDAC1 is subject to both phosphorylation and acetylation (PMID: 23233904). The apparent molecular weight of VDAC1 is 30-37 kDa (PMID: 14573604; 23754752; 25681439). Hypoxic conditions were found to trigger cleavage of the VDAC1 C-terminal to yield a 26-kDa truncated but active form (PMID: 22389449; 23233904). This antibody is specific to VDAC1.

Lagerung

Lagerungsbedingungen:

Bei -20°C lagern. Vor Licht schützen. Nach dem Versand ein Jahr stabil.

Lagerungspuffer:

BS mit 50% Glycerin, 0,05% Proclin300, 0,5% BSA, pH 7,3.

Aliquotieren ist nicht notwendig bei -20°C Lagerung

*** 20ul-Größen enthalten 0.1% BSA


For technical support and original validation data for this product please contact:

T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or 1(312) 455-8498 (outside USA)

E: proteintech@ptglab.com
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Ausgewählte Validierungsdaten

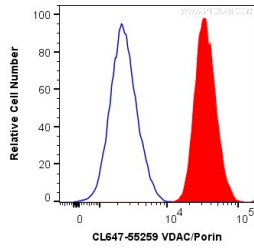

1X10⁶ HepG2 cells were intracellularly stained with 0.2 µg CoraLite® Plus 647 Anti-Human VDAC1/Porin (CL647-55259) (red), or 0.2 µg Isotype Control (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).