

À des fins de recherche uniquement

# Anticorps Polyclonal de lapin anti-Collagen Type VI

Numéro de catalogue: 14853-1-AP

7 Publications



## Informations de base

Numéro de catalogue:	BC065509	Méthode de purification:
14853-1-AP	1292	Purification par affinité contre l'antigène
Taille:	Identification du gène (NCBI):	Dilutions recommandées:
150ul , Concentration: 500 µg/ml by Nanodrop and 267 µg/ml by Bradford method using BSA as the standard;	collagen, type VI, alpha 2	WB 1:1000-1:4000 IHC 1:500-1:2000
Hôte:	MW calculé	
Lapin	109 kDa	
Isotype:	MW observés:	
IgG	140 kDa	
Immunogen Catalog Number:		
AG6635		

## Applications

Applications testées:	Contrôles positifs:
IHC, WB, ELISA	WB : tissu placentaire humain, tissu cardiaque humain, tissu cutané de porc
Demandes citées:	IHC : tissu d'estomac humain, tissu de cancer du foie humain
IF, IHC, WB	
Spécificité de l'espèce:	
Humain, porc	
Espèces citées:	
Humain	
<i>Remarque-IHC: il est suggéré de démasquer l'antigène avec un tampon de TE buffer pH 9,0; (*) À défaut, 'le démasquage de l'antigène peut être effectué avec un tampon citrate pH 6,0.</i>	

## Informations générales

Collagen VI, composed of three polypeptides, alpha-1(VI), alpha-2(VI), and alpha-3(VI) chains, is a beaded filament collagen found in most connective tissues. The alpha-1(VI) and alpha-2(VI) chains are 140 kDa, while the alpha-3(VI) chain is 200-300 kDa (PMID: 3792302, 3348212). COL6A2 gene encodes the alpha-2 subunit of type VI collagen. The alpha-2(VI) subunit has a calculated molecular weight of 109 kDa, and an apparent molecular weight of 140 kDa due to posttranslational modification. Mutations in COL6A2 gene are associated with Bethlem myopathy and Ullrich scleroatonic muscular dystrophy. Three transcript variants have been identified for this gene.

## Publications notables

Autrice	Pubmed ID	Journal	Application
Emna Ouni	35341935	Matrix Biol	IF
Naomi A Carne	30890563	Mol Cell Proteomics	WB
Naomi A. Carne	33526191	Mol Cell Proteomics	WB

## Stockage

### Stockage:

Stocker à -20°C. Stable pendant un an après l'expédition.

### Tampon de stockage:

PBS avec azoture de sodium à 0,02 % et glycérol à 50 % pH 7,3

L'aliquotage n'est pas nécessaire pour le stockage à -20C


\*\*\* Les 20ul contiennent 0,1% de BSA.

For technical support and original validation data for this product please contact:  
T: 1(888) 4PTGLAB (1-888-478-4522) (toll free  
in USA), or 1(312) 455-8498 (outside USA)


E: proteintech@ptglab.com  
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


## Données de validation sélectionnées

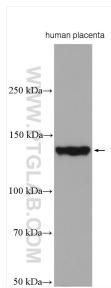

human placenta tissue were subjected to SDS PAGE followed by western blot with 14853-1-AP (Collagen Type VI antibody) at dilution of 1:4000 incubated at room temperature for 1.5 hours.

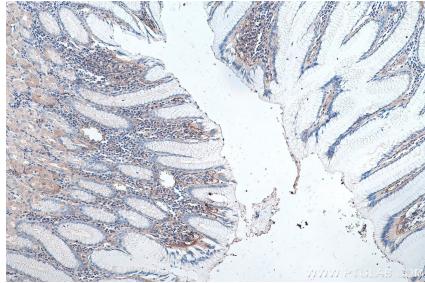

Immunohistochemical analysis of paraffin-embedded human stomach tissue slide using 14853-1-AP (Collagen Type VI antibody) at dilution of 1:1000 (under 10x lens). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0).

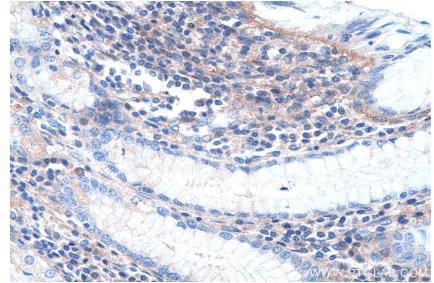

Immunohistochemical analysis of paraffin-embedded human stomach tissue slide using 14853-1-AP (Collagen Type VI antibody) at dilution of 1:1000 (under 40x lens). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0).