

À des fins de recherche uniquement

Anticorps Polyclonal de lapin anti-PRDX4



Numéro de catalogue: **CL488-10703**

Phare

Informations de base

Numéro de catalogue: CL488-10703	Numéro d'acquisition GenBank: BC007107	Méthode de purification: Purification par affinité contre l'antigène
Taille: 100ul, Concentration: 1000 µg/ml by Nanodrop;	Identification du gène (NCBI): 10549	Dilutions recommandées: IF 1:50-1:500
Hôte: Lapin	Nom complet: peroxiredoxin 4	Excitation/Emission maxima wavelengths: 493 nm / 522 nm
Isotype: IgG	MW calculé 31 kDa	
Immunogen Catalog Number: AG1177	MW observés: 30 kDa	

Applications

Applications testées: FC (Intra), IF	Contrôles positifs: IF : cellules HepG2,
Spécificité de l'espèce: Humain, rat, souris	

Informations générales

PRDX4 (Peroxiredoxin-4) is also named as AOE37-2, Prx-IV and belongs to the AhpC/TSA family. PRDX4 is associated with acrosome formation during rat spermatogenesis and has a protective role in the male reproductive tract, because the phenotype of mice lacking this isoform includes testicular atrophy and increased sperm DNA damage (PMID:20864641). It is initially synthesized as a membrane-binding 31-kDa protein and processed into a 27-kDa secretory form and is discarded with the residual bodies (PMID:19208552). PRDX4 is a pentamer of dimers (PMID:23025503). This antibody is specific to PRDX4.

Stockage


Stockage:
Stocker à -20 °C. Éviter toute exposition à la lumière. Stable pendant un an après l'expédition.
Tampon de stockage:
PBS avec glycérol à 50 %, Proclin300 à 0,05 % et BSA à 0,5 %, pH 7,3.
L'aliquotage n'est pas nécessaire pour le stockage à -20C

***** Les 20ul contiennent 0,1% de BSA.**


For technical support and original validation data for this product please contact:
T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or 1(312) 455-8498 (outside USA) E: proteintech@ptglab.com W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Données de validation sélectionnées

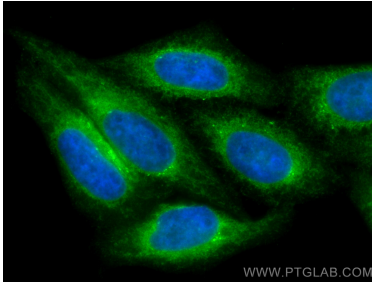

Immunofluorescent analysis of (4% PFA) fixed HepG2 cells using CoraLite® Plus 488 PRDX4 antibody (CL488-10703) at dilution of 1:200.

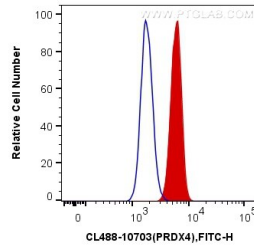

1X10⁶ HepG2 cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human PRDX4 (CL488-10703) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).