

À des fins de recherche uniquement

Anticorps Polyclonal de lapin anti-KAP1

Numéro de catalogue:**CL488-15202** Phare



Informations de base

Numéro de catalogue:	Numéro d'acquisition GenBank:	Méthode de purification:
CL488-15202	BC004978	Purification par affinité contre l'antigène
Taille:	Identification du gène (NCBI):	Excitation/Emission maxima wavelengths:
100ul , Concentration: 1000 µg/ml by Nanodrop;	10155	493 nm / 522 nm
Hôte:	Nom complet:	
Lapin	tripartite motif-containing 28	
Isotype:	MW calculé	
IgG	89 kDa	
Immunogen Catalog Number:	MW observés:	
AG7350	100 kDa	

Applications

Applications testées:

FC (Intra)

Spécificité de l'espèce:

Humain

Informations générales

KAP1, also named as TRIM28 or RNF96, is a 835 amino acid protein, which contain one RING-type zinc finger, one PHD-type zinc finger, one bromo domain and two B box-type zinc fingers. KAP1 localizes in the nucleus and Belongs to the TRIM/RBCC family. KAP1 is a nuclear corepressor for KRAB domain-containing zinc finger proteins and mediates gene silencing by recruiting CHD3, a subunit of the nucleosome remodeling and deacetylation (NuRD) complex, and SETDB1 to the promoter regions of KRAB target genes. KAP1 is expressed in all tissues tested including spleen, thymus, prostate, testis, ovary, small intestine, colon and peripheral blood leukocytes. The calculated molecular weight of KAP1 is 89 kDa, but modified KAP1 is about 100 kDa. (PMID: 18590578)

Stockage

Stockage:

Stocker à -20°C. Éviter toute exposition à la lumière.

Tampon de stockage:

PBS avec glycérol à 50 %, Proclin300 à 0,05 % et BSA à 0,5 %, pH 7,3.


L'aliquotage n'est pas nécessaire pour le stockage à -20°C

*** Les 20ul contiennent 0,1% de BSA.

For technical support and original validation data for this product please contact:
T: 1(888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA) E: proteintech@ptglab.com
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Données de validation sélectionnées

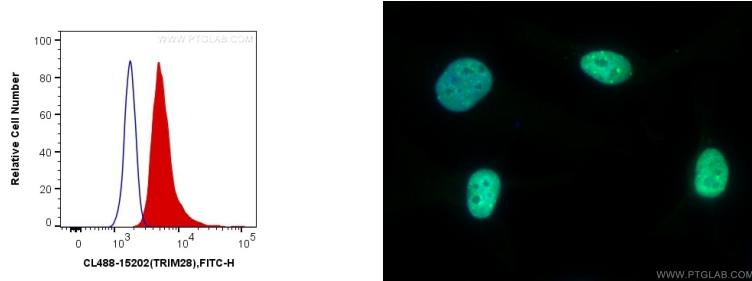

1X10⁶ HeLa cells were intracellularly stained with 0.4 ug Coralite® Plus 488 Anti-Human KAP1 (CL488-15202) (red), or 0.4 ug Control Antibody (CL488-15202) at dilution of 1:200. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).

Immunofluorescent analysis of (4% PFA) fixed HeLa cells using Coralite® Plus 488 KAP1 antibody (CL488-15202) at dilution of 1:200.