

À des fins de recherche uniquement

# Anticorps Monoclonal anti-GORASP2

Numéro de catalogue: CL488-66627



## Informations de base

Numéro de catalogue: CL488-66627	Numéro d'acquisition GenBank: BC007770	Méthode de purification: Purification par protéine A
Taille: 100ul , Concentration: 1000 µg/ml by Nanodrop;	Identification du gène (NCBI): 26003	CloneNo.: 1C9A3
Hôte: Mouse	Nom complet: golgi reassembly stacking protein 2, 55kDa	Dilutions recommandées: IF 1:50-1:500
Isotype: IgG1	MW calculé 55 kDa	Excitation/Emission maxima wavelengths: 493 nm / 522 nm
Immunogen Catalog Number: AG17607	MW observés: 59 kDa	

## Applications

Applications testées: FC (Intra), IF	Contrôles positifs: IF : cellules HepG2,
Spécificité de l'espèce: Humain, porc, rat	

## Informations générales

GORASP2 (Golgi reassembly-stacking protein 2, or GRS2) has multiple phosphorylation sites and can be phosphorylated by MKK/ERK pathway during mitosis. GORASP2 plays a role in the assembly and membrane stacking of the Golgi cisternae, and in the process by which Golgi stacks reform after mitotic breakdown. GORASP2 is supposed to regulate the intracellular transport and presentation of a defined set of transmembrane proteins, such as transmembrane TGFA.

## Stockage

**Stockage:**  
Stocker à -20 °C. Éviter toute exposition à la lumière.  
**Tampon de stockage:**  
PBS avec glycérol à 50 %, Proclin300 à 0,05 % et BSA à 0,5 %, pH 7,3.  
L'aliquotage n'est pas nécessaire pour le stockage à -20C

**\*\*\* Les 20ul contiennent 0,1% de BSA.**


For technical support and original validation data for this product please contact:

T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free  
in USA), or 1(312) 455-8498 (outside USA)


E: [proteintech@ptglab.com](mailto:proteintech@ptglab.com)  
W: [ptglab.com](http://ptglab.com)

**This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.**

## Données de validation sélectionnées

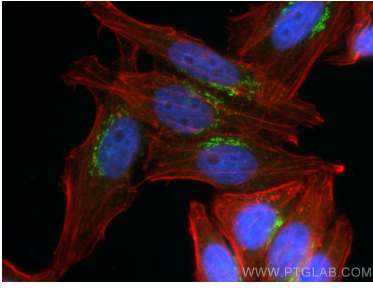

Immunofluorescent analysis of (4% PFA) fixed HepG2 cells using CoraLite® Plus 488 GORASP2 antibody (CL488-66627, Clone: 1C9A3) at dilution of 1:200, CL594-Phalloidin (red).

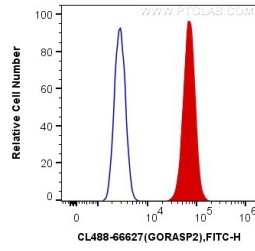

1X10<sup>6</sup> HepG2 cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human GORASP2 (CL488-66627, Clone:1C9A3) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).