

À des fins de recherche uniquement

# Anticorps Polyclonal de lapin anti-Cytokeratin 14

Numéro de catalogue:**CL647-10143** Phare



## Informations de base

Numéro de catalogue:	BC002690	Méthode de purification:
CL647-10143		Purification par affinité contre l'antigène
Taille:	Identification du gène (NCBI):	Excitation/Emission maxima wavelengths:
100ul , Concentration: 1000 µg/ml by Nanodrop;	3861	654 nm / 674 nm
Hôte:	Nom complet:	
Lapin	keratin 14	
Isotype:	MW calculé	
IgG	472 aa, 52 kDa	
Immunogen Catalog Number:	MW observés:	
AG0188	47-50 kDa	

## Applications

### Applications testées:

FC (Intra)

### Spécificité de l'espèce:

Humain, rat, souris

## Informations générales

Keratins are a large family of proteins that form the intermediate filament cytoskeleton of epithelial cells, which are classified into two major sequence types. Type I keratins are a group of acidic intermediate filament proteins, including K9-K23, and the hair keratins Ha1-Ha8. Type II keratins are the basic or neutral counterparts to the acidic type I keratins, including K1-K8, and the hair keratins, Hb1-Hb6. Keratin 14 is a type I cytokeratin. It is usually found as a heterotetramer with keratin 5. Keratins K14 and K5 have long been considered to be biochemical markers of the stratified squamous epithelia, including epidermis.

## Stockage

### Stockage:

Stocker à -20 °C. Éviter toute exposition à la lumière. Stable pendant un an après l'expédition.

### Tampon de stockage:

PBS avec glycérol à 50 %, Proclin300 à 0,05 % et BSA à 0,5 %, pH 7,3.


L'aliquotage n'est pas nécessaire pour le stockage à -20°C

\*\*\* Les 20ul contiennent 0,1% de BSA.

For technical support and original validation data for this product please contact:  
T: 1(888) 4PTGLAB (1-888-478-4522) (toll free  
in USA), or 1(312) 455-8498 (outside USA)      E: proteintech@ptglab.com  
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Données de validation sélectionnées

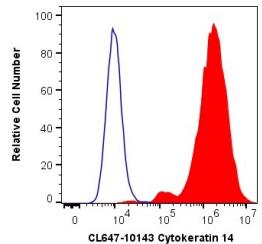

$1 \times 10^6$  A431 cells were intracellularly stained with 0.2 ug Coralite® Plus 647 Anti-Human Cytokeratin 14 (CL647-10143) (red), or 0.2 ug isotype control (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).