For Research Use Only

APOE Monoclonal antibody, PBS Only (Capture)




Catalog Number:68696-1-PBS

Basic Information	Catalog Number: 68696-1-PBS	GenBank Accession Number: BC003557	Purification Method: Protein G purification
	Size: 100ug , Concentration: 1mg/ml by Nanodrop; Source: Mouse Isotype: IgG1 Immunogen Catalog Number: EG0448	GeneID (NCBI): 348 UNIPROT ID: P02649 Full Name: Apolipoprotein E Calculated MW: 36 kDa	CloneNo.: 1A3C5
Applications	Tested Applications: Sandwich ELISA, Indirect ELISA, San Species Specificity: human	nple test	
Product Information	68696-1-PBS targets APOE as part of	f a matched antibody pair:	
	MP50078-1: 68696-1-PBS capture and 68696-2-PBS detection (validated in Sandwich ELISA)		
	Unconjugated mouse monoclonal antibody pair in PBS only (BSA and azide free) storage buffer at a concentration of 1 mg/mL, ready for conjugation.		
	This conjugation ready format makes antibodies ideal for use in many applications including: ELISAs, multiplex assays requiring matched pairs, mass cytometry, and multiplex imaging applications.Antibody use should be optimized by the end user for each application and assay.		
Storage	Storage: Store at -80°C. Storage Buffer: PBS only, pH7.3		


For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA)E: proteintech@ptglab.comW: ptglab.comW: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


Selected Validation Data

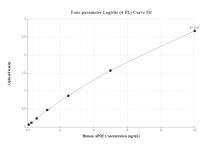


Sandwich ELISA standard curve of MP50078-1, APOE Monoclonal Matched Antibody Pair, PBS Only. Capture antibody: 68696-1-PBS. Detection antibody: 68696-2-PBS. Standard: Eg0448. Range: 0.156-10 ng/mL

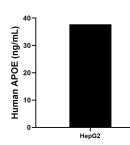


HepG2 cells were cultured in DMEM supplemented with 10% fetal bovine serum, 2 mM L-glutamine, 100 U/mL penicillin, and 100 µg/mL streptomycin sulfate. An aliquot of the cell culture supernate was removed, assayed for human APOE, and measured 37.72 ng/mL

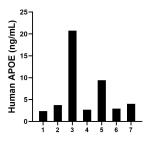

Urine of seven individual healthy human donors was measured. The APOE concentration of detected samples was determined to be 6.56 ng/mL with a range of 2.35-20.75 ng/mL

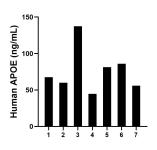

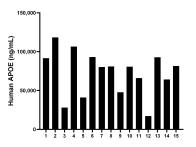

Human milk of seven individual healthy human donors was measured. The APOE concentration of detected samples was determined to be 76.11 ng/mL with a range of 44.74-137.34 ng/mL Serum of sixteen individual healthy human donors was measured. The APOE concentration of detected samples was determined to be 72,675.66 ng/mL with a range of 17,060.26-118,380.99 ng/mL