
FOR IN VITRO RESEARCH USE ONLY. NOT FOR USE IN HUMANS OR ANIMALS.

Recombinant human COQ5 protein



Basic Information	Catalog Number:	Peptide Sequence:
	Ag10203	MAAPGSCALWSYCGRGWSRAMRGCQLLGLRSSWPC
	Size:	
	50 µg	ESVAKKYDVMNDMMSLGIHRVWKDLLLWKMHPLPGT QLLDVAGGTGDIAFRFLNYVQSQHQRKQKRQLRAQ
	Form:	QNLSWEEIAKEYQNEEDSLGGSRVVVCDINKEMLKV
	Available lyophilized	GKQKALAQGYRAGLAWVLGDAEELPFDDDKFDIYTIA
	buman	FGIRNVTHIDQALQEAHRVLKPGGRFLCLEFSQVNNF
		LISRLYDLYSFQVIPVLGEVIAGDWKSYQYLVESIRRFP
		(1 527 00 chooce by be 107074)
	Not tested	
	Endotoxin Level:	
	Please contact the lab for more information	
Reconstitution	Reconstitution:	Shipping:
and Storage	Reconstitute at 0.25 $\mu g/\mu l$ in 200 μl sterile water for short-term	The product is shipped at ambient temperature.
	storage. After reconstitution with sterile water, if glycerol has no effect	Upon receipt, store it immediately at the
	on subsequent experiments, it is recommended to add an	recommended temperature (see below).
	equal volume of glycerol for long-term storage (see Stability	
	and Storage for more details). If a different concentration is needed for your purposes please	
	adjust the reconstitution volume as required (please note: the	
	ion concentration of the final solution will vary according to	
	the volume used). Note: Centrifuge vial before opening. When reconstituting,	
	gently pipet and wash down the sides of the vial to ensure full	
	recovery of the protein into solution.	
Purity	85%, by SDS-PAGE with Coomassie Brilliant Blue staining.	
Formulation	The purified protein was Lyophilized from sterile PBS (58mM Na2HPO4, 17mM NaH2PO4, 68mM NaCl, pH7.4). 5 % trehalose and 5 % mannitol are added as protectant before lyophilization. The elution buffer contain 300mM	
	imidazole.	
Stability and Storage	Store for up to 12 months at -20°C to -80°C as lyophilized powder.	
Storage of	Short Term Storage:	
Reconstituted Protein	Store at 2-8°C for (1-2 weeks). Long Term Storage:	
	Aliquot and store at -20°C to -80°C for up to 3 months,	
	reconstitution with sterile water and addition of an equal volume of glycerol. Avoid repeat freeze-thaw cycles.	
	volume of grycerol. Avoid repeat neeze-ulaw cycles.	
Selected Validation Data		
	74 kBa→	

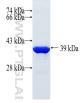

For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522)E : proteintech@ptglab.com W:(toll free in USA), or 1(312) 455-8498ptglab.com(outside USA)(total free in USA)

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.