For Research Use Only

CoraLite® Plus 488-conjugated MOV10 Polyclonal antibody

www.ptglab.com

Purification Method:

IF/ICC 1:50-1:500

wavelengths: 493 nm / 522 nm

Antigen affinity purification

Excitation/Emission maxima

Recommended Dilutions:

Catalog Number: CL488-10370

Featured Product

Basic Information

Catalog Number: GenBank Accession Number:

CL488-10370 BC002548 GeneID (NCBI):

100ul, Concentration: 1000 ug/ml by 4343 Nanodrop:

UNIPROT ID: Q9HCE1 Rabbit Full Name:

Isotype: Mov10, Moloney leukemia virus 10,

IgG homolog (mouse) Immunogen Catalog Number: Calculated MW: AG0384 114 kDa Observed MW:

110-115 kDa

Applications

Tested Applications: IF/ICC, FC (Intra)

Species Specificity: human, mouse, rat

Positive Controls:

IF/ICC: HepG2 cells,

Background Information

MOV10, also named as KIAA1631, belongs to the DNA2/NAM7 helicase family and SDE3 subfamily. It is required for RNA-mediated gene silencing by the RNA-induced silencing complex (RISC). Human MOV10 may regulate a wide range of RNA viruses and could also control the retrotransposition of endogenous retroelements in mammals(PMID:22727223). MOV10 has a broad antiretroviral activity that can target a wide range of retroviruses, and it could be actively involved in host defense against retroviral infection. MOV 10 can potently inhibit HIV-1 replication at multiple stages(PMID:20215113). It is involved in the progression of telomerase-catalyzing reaction via the interaction of telomerase protein and telomere DNA.

Storage


Storage:

Store at -20°C. Avoid exposure to light. Stable for one year after shipment.


PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

Selected Validation Data

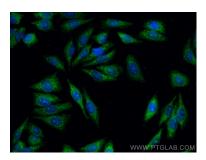

Immunofluorescent analysis of (-20°C Methanol) fixed HepG2 cells using CoraLite® Plus 488 MOV 10 antibody (CL488-10370) at dilution of 1:200.

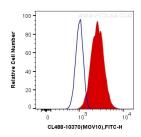

1X10^6 HeLa cells were intracellularly stained with 0.4 ug Coralite® Plus 488 Anti-Human MOV10 (CL488-10370) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).