

For Research Use Only

# CoraLite® Plus 488-conjugated DIABLO Polyclonal antibody



Catalog Number: CL488-10434

Featured Product

## Basic Information

<b>Catalog Number:</b> CL488-10434	<b>GenBank Accession Number:</b> BC011909	<b>Purification Method:</b> Antigen affinity purification
<b>Size:</b> 100ul , Concentration: 1000 µg/ml by Nanodrop;	<b>GeneID (NCBI):</b> 56616	<b>Recommended Dilutions:</b> IF 1:50-1:500
<b>Source:</b> Rabbit	<b>Full Name:</b> diablo homolog (Drosophila)	<b>Excitation/Emission maxima wavelengths:</b> 493 nm / 522 nm
<b>Isotype:</b> IgG	<b>Calculated MW:</b> 27 kDa	
<b>Immunogen Catalog Number:</b> AG0696	<b>Observed MW:</b> 20 kDa	

## Applications

<b>Tested Applications:</b> IF	<b>Positive Controls:</b> IF : HeLa cells,
<b>Species Specificity:</b> human, mouse, rat	

## Background Information

Smac/DIABLO is a mitochondrial protein that promotes apoptosis by neutralizing members of the IAP family. SMAC/DIABLO is ubiquitously expressed in normal tissues, and the expression of SMAC messenger RNA (mRNA) in adult testis is the highest. Smac/DIABLO encodes a 239 amino acid protein (~27 kDa). The first 55 amino acids are a mitochondrial targeting signal peptide that is cleaved to generate mature Smac-DIABLO (~21 kDa). When released from mitochondria, Smac-DIABLO acts as a natural antagonist of inhibitor of apoptosis proteins (IAPs). Smac-DIABLO antagonistic activity is based on its N-terminal tetrapeptide (AVPI) that binds baculoviral IAP repeat (BIR) domains of IAPs, releasing their inhibitory effects on both initiator and effector caspases, thus promoting cell death. (PMID: 32325691, PMID: 32575872, PMID: 25650938)

## Storage

**Storage:**  
Store at -20°C. Avoid exposure to light. Stable for one year after shipment.  
**Storage Buffer:**  
PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.  
Aliquoting is unnecessary for -20°C storage

\*\*\* 20ul sizes contain 0.1% BSA


For technical support and original validation data for this product please contact:

T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free  
in USA), or 1(312) 455-8498 (outside USA)


E: [proteintech@ptglab.com](mailto:proteintech@ptglab.com)  
W: [ptglab.com](http://ptglab.com)

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Selected Validation Data

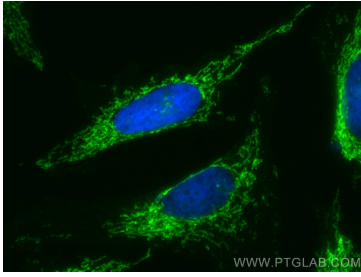

Immunofluorescent analysis of (4% PFA) fixed HeLa cells using CoraLite® Plus 488 DIABLO antibody (CL488-10434) at dilution of 1:100.

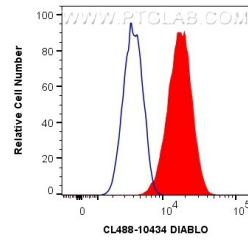

1X10<sup>6</sup> HeLa cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human DIABLO (CL488-10434) (red), or 0.8 ug Isotype Control. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).