For Research Use Only

CoraLite® Plus 488-conjugated Smac/DIABLO Polyclonal antibody Catalog Number:CL488-10434 Featured Product



Basic Information	Catalog Number: CL488-10434	GenBank Accession Number: BC011909	Purification Method: Antigen affinity purification	
	Rabbit	UNIPROT ID: Q9NR28	Recommended Dilutions: IF/ICC 1:50-1:500 Excitation/Emission maxima wavelengths: 493 nm / 522 nm	
				Full Name: diablo homolog (Drosophila) Calculated MW:
		27 kDa		
		Observed MW: 20 kDa		
		Applications	Tested Applications: IF/ICC, FC (Intra)	Positive Controls: IF/ICC : HeLa cells,
Species Specificity: human, mouse, rat				
Background Information	Smac/DIABLO is a mitochondrial protein that promotes apoptosis by neutralizing members of the IAP family. SMAC/DIABLO is ubiquitously expressed in normal tissues, and the expression of SMAC messenger RNA (mRNA) in adult testis is the highest .Smac/DIABLO encodes a 239 amino acid protein (~27 kDa). The first 55 amino acids are a mitochondrial targeting signal peptide that is cleaved to generate mature Smac-DIABLO (~21 kDa). When released from mitochondria, Smac-DIABLO acts as a natural antagonist of inhibitor of apoptosis proteins (IAPs). Smac-DIABLO antagonistic activity is based on its N-terminal tetrapeptide (AVPI) that binds baculoviral IAP repeat (BIR) domains of IAPs, releasing their inhibitory effects on both initiator and effector caspases, thus promoting cell death. (PMID: 32325691, PMID: 32575872, PMID: 25650938)			
Storage	Storage: Store at -20°C. Avoid exposure to light. Stable for one year after shipment. Storage Buffer: PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3. Aliguoting is unnecessary for -20°C storage			


For technical support and original validation data for this product please contact: T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free E: proteintech@ptglab.com in USA), or 1(312) 455-8498 (outside USA) W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data



Immunofluorescent analysis of (4% PFA) fixed HeLa cells using CoraLite® Plus 488 DIABLO antibody (CL488-10434) at dilution of 1:100.

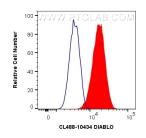

1X10⁶ HeLa cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human DIABLO (CL488-10434) (red), or 0.8 ug Isotype Control. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).