## For Research Use Only

## CoraLite® Plus 488-conjugated PIN1 Polyclonal antibody

Catalog Number:CL488-10495

Featured Product




Basic Information	Catalog Number: CL488-10495	GenBank Accession Number: BC002899	Purification Method: Antigen affinity purification
	Size: 100ul , Concentration: 1000 ug/ml by	GenelD (NCBI): 5300	Recommended Dilutions: IF/ICC 1:50-1:500
	Nanodrop; Source: Rabbit Isotype: IgG Immunogen Catalog Number: AG0767	UNIPROT ID: Q13526 Full Name: peptidylprolyl cis/trans isomerase, NIMA-interacting 1 Calculated MW: 18 kDa Observed MW: 18 kDa	Excitation/Emission maxima wavelengths: 493 nm / 522 nm
Applications	Tested Applications: IF/ICC, FC (Intra) Species Specificity:	Positive Controls: IF/ICC : HEK-293 cells,	
Background Informatior	<ul> <li>PIN1(Peptidyl-prolyl cis-trans isomerase NIMA-interacting 1) is essential for mitosis progression in yeast cells and is hypothesized to perform the same role in mammalian cells. It might regulate cellular processes distinct from the cell cycle itself, such as terminal differentiation through a modulation of differentiation-specific gene expression(PMID:20801874). It colocalizes with NEK6 in the nucleus. Pin1 inhibition simultaneously blocks multiple cancer pathways, disrupts the desmoplastic and immunosuppressive TME, and upregulates PD-L1 and ENT1, rendering pancreatic ductal adenocarcinoma (PDAC) eradicable by immunochemotherapy (PMID: 34388391).</li> </ul>		
Storage	Storage: Store at -20°C. Avoid exposure to ligh Storage Buffer: PBS with 50% Glycerol, 0.05% Procli Aliquoting is unnecessary for -20°C s	n300, 0.5% BSA, pH 7.3.	


For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free<br/>in USA), or 1(312) 455-8498 (outside USA)E: proteintech@ptglab.comW: ptglab.comW: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Selected Validation Data

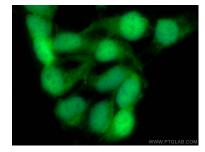

Immunofluorescent analysis of (4% PFA) fixed HEK-293 cells using CoraLite® Plus 488 PIN1 antibody (CL488-10495) at dilution of 1:200.

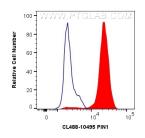

1X10<sup>6</sup> HeLa cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human PlN1 (CL488-10495) (red), or 0.8 ug Isotype Control. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).