#### For Research Use Only

# CoraLite® Plus 488-conjugated RABL3 Polyclonal antibody

Catalog Number:CL488-11916 Featured Product

Antibodies | ELISA kits | Proteins www.ptglab.com

Catalog Number: GenBank Accession Number: **Basic Information** CL488-11916 BC020832 GenelD (NCBI): Size: 100ul , Concentration: 1000 ug/ml by 285282 Nanodrop: UNIPROT ID: Source: Q5HYI8 Rabbit Full Name: Isotype: RAB, member of RAS oncogene lgG family-like 3 Immunogen Catalog Number: Calculated MW: 236 aa, 26 kDa AG2531 **Observed MW:** 26 kDa

#### Purification Method: Antigen affinity purification Recommended Dilutions: IF/ICC 1:50-1:500 Excitation/Emission maxima wavelengths: 493 nm / 522 nm

## **Applications**

Tested Applications: IF/ICC, FC (Intra) Species Specificity: human, mouse Positive Controls:

IF/ICC : A549 cells,

**Background Information** 

The gene of Rabl3 is located on 3q13.3 and is composed of eight exons and seven introns. It is a member of the Rab subfamily, the largest group in the Ras superfamily of small guanosine triphosphatases (GTPases).

**Storage** 

#### Storage:

Store at -20°C. Avoid exposure to light. Stable for one year after shipment. Storage Buffer: PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3. Aliquoting is unnecessary for -20°C storage


For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll freeE: proteintech@ptglab.comin USA), or 1(312) 455-8498 (outside USA)W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

### Selected Validation Data



Immunofluorescent analysis of (-20°C Methanol) fixed A549 cells using CoraLite® Plus 488 RABL3 antibody (CL488-11916) at dilution of 1:200.

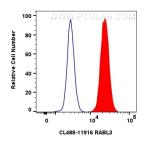

1x10^6 HEK-293 cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human RABL3 (CL488-11916) (red), or 0.8 ug CoraLite® Plus 488conjugated Rabbit IgG control Rabbit PolyAb (CL488-30000, Clone: ) (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).