## CoraLite® Plus 488-conjugated TDRD3 Polyclonal antibody

Catalog Number:CL488-13359 Featured Product

Basic Information	Catalog Number: CL488-13359	GenBank Accession Number: BC030514	Purification Method: Antigen affinity purification
	Size: 100ul , Concentration: 1000 ug/ml by	GenelD (NCBI): 81550	Recommended Dilutions: IF/ICC 1:50-1:500
	Nanodrop; Source: Rabbit Isotype: IgG Immunogen Catalog Number: AG4174	UNIPROT ID: Q9H7E2 Full Name: tudor domain containing 3 Calculated MW: 651 aa, 73 kDa Observed MW: 73 kDa, 83 kDa	Excitation/Emission maxima wavelengths: 493 nm / 522 nm
Applications	Tested Applications: IF/ICC, FC (Intra) Species Specificity: human, mouse, rat	Positive Co IF/ICC : He	
Background Information	TDRD3 (Tudor domain-containing protein 3) contains tudor domain which is an approximately 60-amino acid structure motif. TDRD3 functions as a scaffolding protein that specifically recognizes and binds dimethylarginine-containing proteins. TDRD3 is a transcriptional coactivator that promotes transcription by binding methylarginine marks on histone tails. TDRD3 also possesses an oligosaccharide/nucleotide binding fold and an ubiquitin associated domain capable of binding tetra-ubiquitin, and it is reported to associate with polyribosomes in HeLa cells.		
Storage	Storage: Store at -20°C. Avoid exposure to ligh Storage Buffer: PBS with 50% Glycerol, 0.05% Proclin Aliquoting is unnecessary for -20°C s	n300, 0.5% BSA, pH 7.3.	ıt.

For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll freeE: proteintech@ptglab.comin USA), or 1(312) 455-8498 (outside USA)W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


proteintech

Antibodies | ELISA kits | Proteins WWW.ptglab.com

## Selected Validation Data



Immunofluorescent analysis of (4% PFA) fixed HeLa cells using CoraLite® Plus 488 TDRD3 antibody (CL488-13359) at dilution of 1:200, CL594-Phalloidin (red).

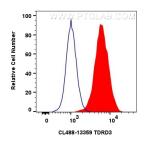

1x10^6 HEK-293 cells were intracellularly stained with 0.8 ug Coralite® Plus 488 Tdrd3 Polyclonal Antibody (CL488-13359)(red), or 0.8 ug Coralite® Plus 488-conjugated Rabbit IgG control Rabbit PolyAb (CL488-30000) (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).