

For Research Use Only

# CoraLite® Plus 488-conjugated PRDX6 Polyclonal antibody

Catalog Number: CL488-13585

Featured Product



## Basic Information

Catalog Number:

CL488-13585

Size:

100ul , Concentration: 1000 ug/ml by Nanodrop;

Source:

Rabbit

Isotype:

IgG

Immunogen Catalog Number:

AG4512

GenBank Accession Number:

BC035857

GeneID (NCBI):

9588

UNIPROT ID:

P30041

Full Name:

peroxiredoxin 6

Calculated MW:

224 aa, 25 kDa

Observed MW:

25-30 kDa

Purification Method:

Antigen affinity purification

Recommended Dilutions:

IF/ICC 1:250-1:1000

Excitation/Emission maxima wavelengths:

493 nm / 522 nm

## Applications

Tested Applications:

IF/ICC, FC (Intra)

Species Specificity:

human, mouse, rat

Positive Controls:

IF/ICC : HeLa cells,

## Background Information

PRDX6 (Peroxiredoxin-6), also named as AOP2 or KIAA0106, is a unique member of the peroxiredoxin family of antioxidants. PRDX6 is highly expressed in liver and protects cells from oxidative damage by reducing H<sub>2</sub>O<sub>2</sub> and various lipid Peroxides (PMID: 17382207). It can form a dimer (PMID: 20500660). PRDX6 is expressed in all major organs, with a particularly high level in lung (PMID: 15890616). Prdx6 is detected at approximately 24 to 28 kDa, and can be monosumoylated with the molecular mass of about 40 kDa (PMID: 24910119).

## Storage

Storage:

Store at -20°C. Avoid exposure to light. Stable for one year after shipment.

Storage Buffer:

PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage


For technical support and original validation data for this product please contact:

T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or 1(312) 455-8498 (outside USA)


E: [proteintech@ptglab.com](mailto:proteintech@ptglab.com)  
W: [ptglab.com](http://ptglab.com)

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Selected Validation Data

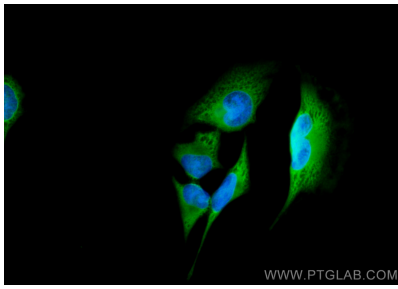

Immunofluorescent analysis of (-20°C Methanol) fixed HeLa cells using CoraLite® Plus 488 PRDX6 antibody (CL488-13585) at dilution of 1:500.

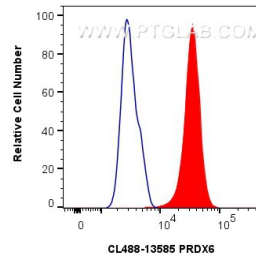

1X10<sup>6</sup> HeLa cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human PRDX6 (CL488-13585) (red), or 0.8 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).