For Research Use Only

CoraLite® Plus 488-conjugated Cytokeratin 19 Polyclonal antibody Catalog Number:CL488-14965 Featured Product


oroteintech Antibodies | ELISA kits | Proteins www.ptglab.com

Basic Information	Catalog Number: CL488-14965	GenBank Accession Number: BC002539	Purification Method: Antigen affinity purification
	Size: 100ul , Concentration: 1000 ug/ml by Nanodrop; Source: Rabbit Isotype: IgG Immunogen Catalog Number: AG6830	GeneID (NCBI): 7 3880 UNIPROT ID: P08727 Full Name: keratin 19 Calculated MW: 44 kDa	Recommended Dilutions: IF/ICC 1:50-1:500 Excitation/Emission maxima wavelengths: 493 nm / 522 nm
Background Information	Keratins are a large family of proteins that form the intermediate filament cytoskeleton of epithelial cells, which are classified into two major sequence types. Type I keratins are a group of acidic intermediate filament proteins, including K9-K23, and the hair keratins Ha1-Ha8. Type II keratins are the basic or neutral courterparts to the acidic type I keratins, including K1-K8, and the hair keratins, Hb1-Hb6. KRT19, one of type I keratins, is specifically expressed in the periderm, the transiently superficial layer that envelopes the developing epidermis. Due to its high sensitivity, KRT19 is the most used marker for detection of tumor cells disseminated in lymph nodes, peripheral blood, and bone marrow of breast cancer patients. This antibody, generated against full length KRT19 protein, also recognizes KRT17, another type I keratin homologous to KRT19.		
Storage	Storage: Store at -20°C. Avoid exposure to ligh Storage Buffer: PBS with 50% Glycerol, 0.05% Proclin Aliquoting is unnecessary for -20°C st	1300, 0.5% BSA, pH 7.3.	nt.


For technical support and original validation data for this product please contact: T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free E: proteintech@ptglab.com in USA), or 1(312) 455-8498 (outside USA) W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data

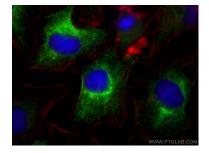

Immunofluorescent analysis of (-20°C Ethanol) fixed HeLa cells using CoraLite® Plus 488 Cytokeratin 19 antibody (CL488-14965) at dilution of 1:200, CL594-phalloidin (red).

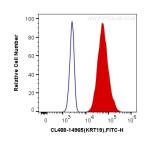

1X10^{^6} MCF-7 cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human Cytokeratin 19 (CL488-14965) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).