

For Research Use Only

# CoraLite®650-conjugated Caspase 9/p35/p10 Polyclonal antibody



Catalog Number:CL650-10380

## Basic Information

<b>Catalog Number:</b> CL650-10380	<b>GenBank Accession Number:</b> BC002452	<b>Purification Method:</b> Antigen affinity purification
<b>Size:</b> 100ul , Concentration: 1000 µg/ml by Nanodrop;	<b>GeneID (NCBI):</b> 842	<b>Excitation/Emission maxima wavelengths:</b> 654 nm / 674 nm
<b>Source:</b> Rabbit	<b>Full Name:</b> caspase 9, apoptosis-related cysteine peptidase	
<b>Isotype:</b> IgG	<b>Calculated MW:</b> 46 kDa	
<b>Immunogen Catalog Number:</b> AG0404	<b>Observed MW:</b> 46 kDa, 35 kDa	

## Applications

**Tested Applications:**  
FC (Intra)

**Species Specificity:**  
human, mouse, rat

## Background Information

Caspase 9, apoptosis-related cysteine protease (CASP9, synonyms: MCH6, APAF3, APAF-3, ICE-LAP6, CASPASE-9c) is a member of the cysteine-aspartic acid protease (caspase) family. Sequential activation of caspases plays a central role in the execution-phase of cell apoptosis. Caspases exist as inactive proenzymes which undergo proteolytic processing at conserved aspartic residues to produce 2 subunits, large and small, that dimerize to form the active enzyme. Caspase 9 is processed by APAF1; this step is thought to be one of the earliest in the caspase activation cascade. 10380-1-AP can recognize the pre- and cleaved- caspase 9.

## Storage

**Storage:**  
Store at -20°C. Avoid exposure to light.


**Storage Buffer:**  
PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

For technical support and original validation data for this product please contact:  
T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or 1(312) 455-8498 (outside USA)      E: proteintech@ptglab.com  
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Selected Validation Data

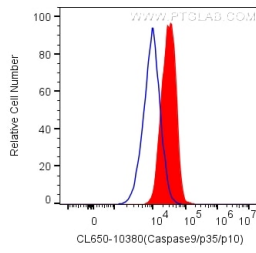

1X10<sup>6</sup> HepG2 cells were intracellularly stained with 0.2 ug CoraLite®650 Anti-Human Caspase 9/p35/p10 (CL650-10380) (red), or 0.2 ug Isotype Control (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).