

## Human ICAM3/CD50 Sandwich ELISA Kit Datasheet

Please read it entirely before use

**Catalogue Number:** KE00633

**Size:** 96T

**Sensitivity:** 0.02 ng/mL

**Range:** 0.313-20 ng/mL

**Usage:** For the quantitative detection of human ICAM3/CD50 concentrations in serum and plasma.


This product is for research use only and not for use in human or animal therapeutic or diagnostic.

<b>Table of content</b>	<b>page</b>
1. Background .....	3
2. Principle .....	3
3. Required Materials .....	3
4. Kit Components and Storage .....	4
5. Safety Notes .....	4
6. Sample Collection and Storage .....	4
7. Regent Preparation .....	5
8. Assay Procedure Summary .....	6
9. Validation Data .....	7
9.1 Standard curve .....	7
9.2 Precision .....	7
9.3 Recovery .....	7
9.4 Sample values .....	8
9.5 Sensitivity .....	8
9.6 Linearity .....	8
9.7 Specificity .....	8
10. References .....	8

## 1. Background

ICAM3 (Intercellular adhesion molecule 3) is also named as CD50, CDw50 and ICAM-R. ICAM3 is a member of the intercellular adhesion molecule (ICAM) family. ICAM3, being expressed constitutively on T cells and other leukocytes, is a potent signalling molecule and a major ligand in the initiation of T cell-mediated immune responses. ICAM3 is the natural ligand for both DC-SIGN and L-SIGN. ICAM3 is abundantly and constitutively expressed in all leukocytes and is probably the most important ligand for LFA-1 in initiating immune responses, and ICAM3 is a typical type I transmembrane glycoprotein. The promoter region of ICAM3 can be bound by co-transcription factors such as ETs and C/EBP, promoting macrophage maturation as well as monocyte extravasation. ICAM3 is widely expressed not only in lymphocytes such as T and B cells but also in many endothelial cells and dendritic cells.

## 2. Principle

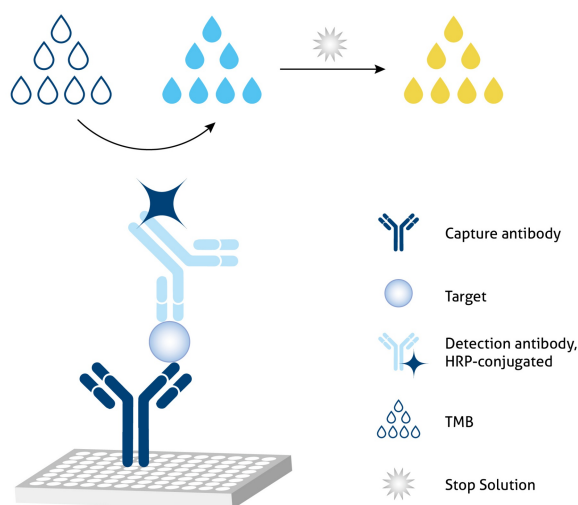

### Sandwich ELISA structure (Detection antibody labeled with HRP)

A capture antibody is pre-coated onto the bottom of wells which binds to analyte of interest. A detection antibody labeled with HRP also binds to the analyte. TMB acts as the HRP substrate and the solution color will change from colorless to blue. A stop solution containing sulfuric acid turns solution yellow. The color intensity is proportional to the quantity of bound protein which is measurable at 450 nm with the correction wavelength set at 630 nm.

## 3. Required Materials

- 3.1 A microplate reader capable of measuring absorbance at 450 nm with the correction wavelength set at 630 nm.
- 3.2 Calibrated, adjustable precision pipettes and disposable plastic tips. A manifold multi-channel pipette is recommended for large assays.
- 3.3 Plate washer: automated or manual.
- 3.4 Absorbent paper towels.
- 3.5 Glass or plastic tubes to prepare standard and sample dilutions.
- 3.6 Beakers and graduated cylinders.
- 3.7 Log-log or semi-log graph paper or computer and software for ELISA data analysis. A four-parameter logistic (4-PL) curve-fit is recommended.

## 4. Kit Components and Storage

Microplate - antibody coated 96-well microplate (8 well × 12 strips)	1 plate	<b>Unopened Kit:</b> Store at 2-8°C for 6 months or -20°C for 12 months. <b>Opened Kit:</b> All reagents stored at 2-8°C for 7 days. <b>Please use a new standard for each assay.</b>
Protein standard - 40 ng/bottle; lyophilized	2 bottles	
Detection antibody, HRP-conjugated (100×) - 120 µL/vial*	1 vial	
Sample Diluent PT 4B1 - 30 mL/bottle	1 bottle	
Detection Diluent - 30 mL/bottle	1 bottle	
Wash Buffer Concentrate (20×) - 30 mL/bottle	1 bottle	
Tetramethylbenzidine Substrate (TMB) - 12 mL/bottle	1 bottle	
Stop Solution - 12 mL/bottle	1 bottle	
Plate Cover Seals	4 pieces	

\* Centrifugation immediately before use

## 5. Safety Notes

- 5.1 Avoid any skin and eye contact with Stop Solution and TMB. In case of contact, wash thoroughly with water.
- 5.2 Do not use the kit after the expiration date.
- 5.3 Do not mix or substitute reagents or materials from other kit lots or other sources.
- 5.4 Be sure to wear protective equipment such as gloves, masks and goggles during the experiment.
- 5.5 When using an automated plate washer, adding a 30 second soak period following the addition of Wash Buffer to improve assay precision

## 6. Sample Collection and Storage

- 6.1 Serum: Allow blood samples to clot for 30 minutes, followed by centrifugation for 15 minutes at 1000xg. Clear serum can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.
- 6.2 Plasma: Use EDTA, heparin, or citrate as an anticoagulant for plasma collection. Centrifuge for 15 minutes at 1000xg within 30 minutes of collection. The plasma can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.

## 7. Regent Preparation

**7.1 Wash Buffer (1X):** If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Add 30 mL of Wash Buffer Concentrate(20X) to 570 mL deionized or distilled water to prepare 1X Wash Buffer.


**7.2 Detection Antibody, HRP-conjugated(1X):** Dilute 100X Detection Antibody, HRP-conjugated 1:100 using Detection Diluent prior to assay. Suggested 1:100 dilution: 10  $\mu$ L 100X Detection Antibody, HRP-conjugated + 990  $\mu$ L Detection Diluent (Centrifuge the 100X Detection Antibody solution, HRP-conjugated for a few seconds prior to use).

**7.3 Sample Dilution:** Different samples should be diluted with corresponding Sample Diluent, samples may require further dilution if the readout values are higher than the highest standard OD reading. Variations in sample collection, processing and storage may affect the results of the measurement.

Recommended Dilution for different sample types: 1:40 or 1:80 is recommended for human serum and plasma.

### 7.4 Standard Serial Dilution:

Add 2 mL Sample Diluent PT 4B1 in protein standard.

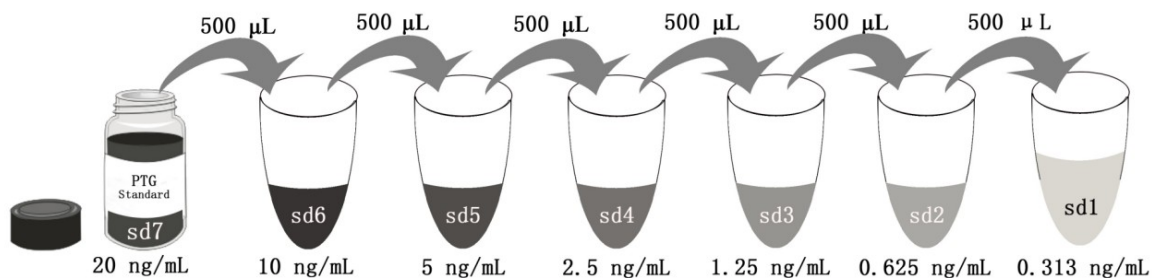

Add # $\mu$ L of Standard diluted in the previous step	—	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L
# $\mu$ L of Sample Diluent PT 4B1	<b>2000 <math>\mu</math>L</b>	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L	500 $\mu$ L
	"sd7"	"sd6"	"sd5"	"sd4"	"sd3"	"sd2"	"sd1"

## 8. Assay Procedure Summary

Bring all reagents to room temperature before use (Detection antibody, HRP-conjugated can be used immediately). To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

8.1 Take out the required number of microplate strips and return excess strips to the foil pouch containing the drying reagent pack and reseal; store at 4°C immediately. Microplate strips should be used in one week.

8.2 Preset the layout of the microplate, including control group, standard group and sample group, add 100 µL of each standard and sample to the appropriate wells. (Make sure sample addition is uninterrupted and completed within 5 to 10 minutes, it is recommended to assay all standards, controls, and samples in duplicate).

8.3 Seal plate with cover seal, pressing it firmly onto top of microwells. Incubate the plate for 2 hours at 37°C.

### 8.4 Wash

1) Gently remove the cover seal. Discard the liquid from wells by aspirating or decanting. Remove any residual solution by tapping the plate a few times on fresh paper towels.

2) Wash 4 times with 1X Wash Buffer, using at least 350-400 µL per well. Following the last wash, firmly tap plates on fresh towels 10 times to remove residual Wash Buffer. Avoid getting any towel fibers in the wells or wells drying out completely.

8.5 Add 100 µL of 1X Detection antibody, HRP-conjugated solution (refer to Reagent Preparation 7.2) to each well. Seal plate with cover seal and incubate for 40 minutes at 37°C.

8.6 Repeat wash step in 8.4.

8.7 Signal development: Add 100 µL of TMB substrate solution to each well, protected from light. Incubate for 15 to 20 minutes. Substrate Solution should remain colorless until added to the plate.

8.8 Quenching color development: Add 100 µL of Stop Solution to each well in the same order as addition of the TMB substrate. Mix by tapping the side of the plate gently. NB: Avoid skin and eye contact with the Stop solution.

8.9 Read results: Immediately after adding Stop solution read the absorbance on a microplate reader at a wavelength of 450 nm. If possible, perform a double wavelength readout (450 nm and 630 nm).


8.10 Data analysis: Calculate the average of the duplicate readings (OD value) for each standard and sample, and subtract the average of the zero standard absorbance. Construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis, use four-parameter logistic curve-fit (4-PL) analysis to do this. If the samples have been diluted, the OD readout from the standard curve must be multiplied by the dilution factor used.

Step	Reagent	Volume	Incubation	Wash	Notes
1	Standard and Samples	100 µL	120 min	4 times	Cover Wells incubate at 37°C
2	Diluent Detection antibody, HRP-conjugated Solution	100 µL	40 min	4 times	Cover Wells incubate at 37°C
3	TMB Substrate	100 µL	15-20 min	Do not wash	Incubate in the dark at 37°C
4	Stop Solution	100 µL	0 min	Do not wash	-
5	Read plate at 450 nm and 630 nm immediately after adding Stop solution. DO NOT exceed 5 minutes.				

## 9. Validation Data

### 9.1 Standard curve

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

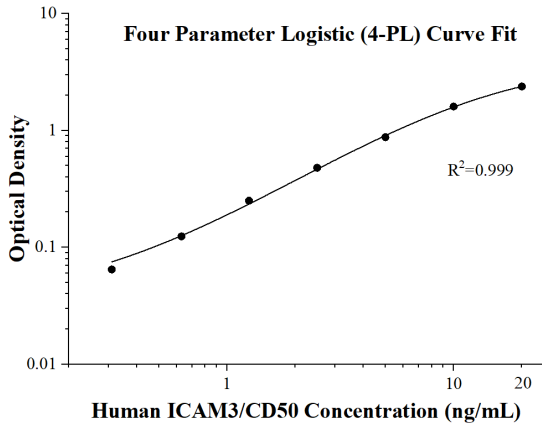

(ng/mL)	O.D	Average	Corrected
0	0.0255 0.0279	0.0267	-
0.313	0.089 0.0938	0.0914	0.0647
0.625	0.1513 0.1503	0.1508	0.1241
1.25	0.2791 0.2756	0.27735	0.25065
2.5	0.492 0.5201	0.50605	0.47935
5	0.8429 0.9616	0.90225	0.87555
10	1.6015 1.6618	1.63165	1.60495
20	2.3636 2.4474	2.4055	2.3788

### 9.2 Precision

**Intra-assay Precision** (Precision within an assay) Three samples of known concentration were tested 8 times on one plate to assess intra-assay precision.

**Inter-assay Precision** (Precision between assays) Three samples of known concentration were tested in 16 separate assays to assess inter-assay precision.

Intra-assay Precision				
Sample	n	Mean (ng/mL)	SD	CV%
1	8	9.72	0.18	1.88
2	8	2.44	0.03	1.36
3	8	1.28	0.04	2.81

Inter-assay Precision				
Sample	n	Mean (ng/mL)	SD	CV%
1	16	9.42	0.41	4.40
2	16	2.32	0.13	5.80
3	16	1.26	0.05	3.89

### 9.3 Recovery

The recovery of human ICAM3/CD50 spiked to three different levels throughout the range of the assay in human serum was evaluated.

Sample Type		Average% of Expected	Range (%)
Human serum	1:160	85	82-89
	1:320	87	82-90

## 9.4 Sample values

**Human serum** - human serum samples were evaluated for the presence of human ICAM3/CD50 in this assay.

Sample Type	Mean (ng/mL)	Range (ng/mL)
Human serum (n=16)	56.07	42.77-80.50

## 9.5 Sensitivity

The minimum detectable dose of human ICAM3/CD50 is 0.02 ng/mL. This was determined by adding two standard deviations to the concentration corresponding to the mean O.D. of 20 zero standard replicates.

## 9.6 Linearity

To assess the linearity of the assay, samples were diluted with the appropriate **Sample Diluent** to produce samples with values within the dynamic range of the assay.

(The human serum samples were initially diluted 1:20.)

		Human serum
1:2	Average% of Expected	100
	Range (%)	-
1:4	Average% of Expected	99
	Range (%)	96-102
1:8	Average% of Expected	103
	Range (%)	101-105
1:16	Average% of Expected	96
	Range (%)	94-98

## 9.7 Specificity

This kit specifically recognizes native and recombinant human ICAM3/CD50.

## 10. References

1. Fawcett J, et al. (1992) Nature. 360(6403):481-4.
2. Montoya MC, et al. (2002) Nat Immunol. 3(2):159-68.
3. Chan KY, et al. (2007) J Infect Dis. 196(2):271-80.
4. Shen X, et al. (2025) J Cancer. 16(5):1563-1574.
5. Estecha A, et al. (2012) PLoS One. 7(3):e33313.
6. Shen X, et al. (2025) J Cancer. 16(5):1563-1574.