

For Research Use Only

Spot-Trap® Agarose, Kit for Immunoprecipitation



www.ptglab.com

Catalog Number: etak

Basic Information

Catalog Number:
etak

Applications:
IP, CoIP, ChIP, RIP

Conjugate:
Agarose beads; bead size: ~ 90 µm (cross-linked 4 % agarose beads)

Host:
Alpaca

Type:
Nanobody

Class:
Recombinant

Description

Spot-Trap® Agarose is an affinity bead for immunoprecipitation (IP) of Spot-tagged proteins. It comprises an anti-Spot-Tag® VHH/ Nanobody coupled to agarose beads.

Binding capacity

17.5 µg of recombinant Spot-tagged protein (~30 kDa) per 25 µL bead slurry

Specificity/Target

Spot-Tag® sequence PDRVRAVSHWSS

Elution buffer

SDS sample buffer
10 mM NaOH pH 12

Affinity (K_D)

Dissociation constant K_D of 6 nM

Storage


Storage:
Shipped at ambient temperature. Upon receipt store at 4°C. Stable for one year. Do not freeze!
Storage Buffer:
20% ethanol

For technical support and original validation data for this product please contact


T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or E: proteintech@ptglab.com
1(312) 455-8498 (outside USA) W: www.ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


Selected Validation Data

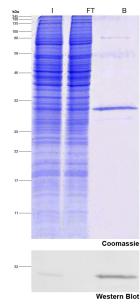

Immunoprecipitation of Spot-tagged protein.

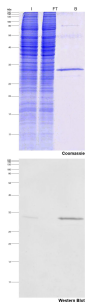

Immunoprecipitation of Spot-tagged protein. Input (I), non-bound (FT) and bound (B), Coomassie staining and Western blotting. Western blot indicates high effectivity.

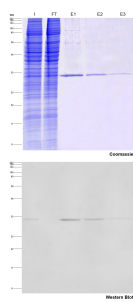

Elution of Spot-tag fusion protein. Input (I), non-bound (FT) and eluted (E1, E2 and E3) fractions. Elutions with Spot-peptide.

	Spot-Trap	Spot-Cap	Spot-Label	Spot V _H H
Immunoprecipitation (IP) / Co-IP / Co-IP/MS	✓			
Protein purification		✓		
Western Blotting			✓	✓
Immunofluorescence			✓	✓

Application of Spot-tag related products.