

For Research Use Only

anti-Rabbit IgG Magnetic VHH Agarose for Immunoprecipitation



www.ptglab.com

Catalog Number: rIGma

Basic Information

Catalog Number:

rIGma

Applications:

IP, Co-IP

Conjugate:

Magnetic Agarose beads; ~40 um (cross-linked 6% magnetic agarose beads)

Host:

Alpaca

Type:

Nanobody

Class:

Recombinant - Animal free production

Description

anti-Rabbit IgG IP Beads is an affinity resin for IP of Rabbit IgG. It consists of rabbit specific VHHs (Nanobodies) coupled to Magnetic Agarose beads.

Binding capacity

/

Elution buffer

SDS Sample Buffer

Wash buffer compatibility

/

Affinity (K_D)

/

Storage

Storage:

+4°C / do not freeze!

Storage Buffer:


20% Ethanol

For technical support and original validation data for this product please contact


T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or E: proteintech@ptglab.com
1(312) 455-8498 (outside USA) W: www.ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data

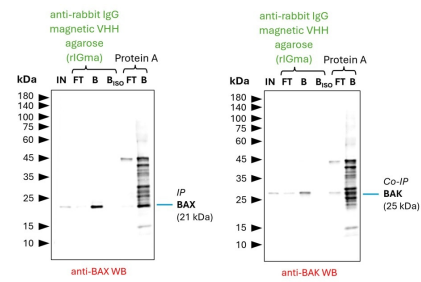

Co-IP of BAX and BAK by anti-rabbit IgG VHH agarose (rIGma) using a BAX polyclonal IgG (Proteintech: 50599-2-Ig). As control an IgG isotype control antibody (Proteintech: 98136-1-RR) was used (BISO). 5 µg of each IgG was spiked into HEK293T cell lysate derived from 0.5×10^7 cells. 1% of input (IN) and flow through (FT) and 25% of bound (B) fraction was loaded onto an SDS-PAGE gel. For Western blot analysis BAX was detected using a polyclonal rabbit IgG (Proteintech: 50599-2-Ig) (1:2000) labeled using FlexAble HRP (Proteintech: KFA005). The presence of BAK co-precipitated with BAX was confirmed using a polyclonal rabbit IgG (Proteintech: 29552-1-AP) (1:2000). While rIGa allows for clean detection of both proteins, Protein A magnetic beads show additional bands, probably due to leaching of Protein A from the beads into the IP fractions, which can lead to binding of the detection antibody and unspecific signals.

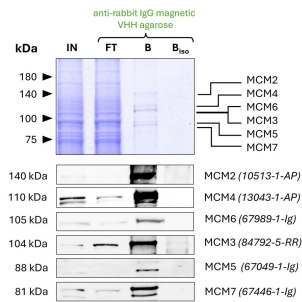

Co-IP of MCM complex via pulldown of MCM6 using 5 µg of polyclonal rabbit anti-MCM6 antibody (PTG: 13347-2-AP) and anti-rabbit IgG magnetic VHH agarose (rIGma). All subunits of the 600 kDa hetero-hexameric complex are successfully precipitated, as shown in a Coomassie blue stained SDS-gel and by Western blot analysis using subunit specific antibodies. Apparent molecular weights are provided. For input (IN) and flowthrough (FT) fractions, 1% was loaded, respectively. For bound (B) fraction and bound fraction of isotype control antibody (BISO PTG: 98136-1-RR), 20% was loaded. Product codes for Proteintech antibodies are provided in parentheses.