For Research Use Only

## Nano-Secondary® anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 568 [CTK0105, CTK0106]



www.ptglab.com

Catalog Number: sms2bAF568-1 1 Publications

Catalog Number: sms2bAF568-1 **Basic Information** 

**Applications:** IF, WB Host: Alpaca Conjugate: Alexa Fluor® 568

Class: Recombinant RRID: AB\_2827582

**Purification Method:** 

Recombinant expression, affinity purification

Type:
Mixture of 2 monoclonal Nanobodies

**Description** 

Nano-Secondary® anti-mouse IgG2b, Fc-specific recombinant VHH is an anti-mouse IgG subclass specific antibody. This secondary antibody product consists of Nanobodies that bind to mouse IgG2b with high affinity & specificity.

**Species Reactivity** 

Mouse IgG2b Fc-fragment No cross-reactivity: goat, guinea pig, rabbit, rat, sheep, human, macaque serum proteins, mouse IgG1, IgG2a, IgG2c, IgG3

**Physical State** 

**Suggested Dilution** 

1(312) 455-8498 (outside USA)

Immunofluorescence 1:1,000

Western blot 1:1,000

Affinity (K<sub>D</sub>)


CTK0105:  $K_D = 5 \text{ nM}$ , CTK0106:  $K_D = 0.2 \text{ nM}$ 

Storage

Storage: Store at -20°C short term or -80°C long term. Aliquot upon delivery. Avoid freeze-thaw cycles.


10 mM HEPES pH 7.0, 500 mM NaCl, 5 mM EDTA, Preservative: 0.09 % Sodium azide

## **Selected Validation Data**

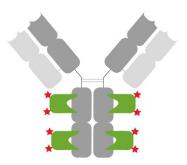

Anti-mouse IgG2b Nano-Secondary: Well-defined and characterized immunostaining. Primary antimouse IgG2b antibody (grey) with 2X2 monoclonal mouse Fc- specific Nanobodies (green) bound. In total, 8 fluorophores (red stars) label the mouse IgG2b primary antibody.



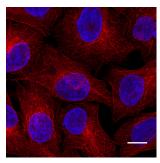

HeLa cells were immunostained with mouse IgG2b anti-\(\beta\)Tubulin antibody + alpaca anti-mouse IgG2b VHH Alexa Fluor\(\text{0}\) 568 (red). Nuclei were stained with DAPI (blue). Scale bar, 10 \(\text{\pm}\)m. Images were recorded at the Core Facility Bioimaging at the Biomedical Center, LMU Munich



One-step immunostaining is the simultaneous incubation of mouse IgG2b primary antibody and anti-mouse IgG2b Nano-Secondary. This method reduces incubation and hands-on time. Simultaneous incubation also supports multiplexing, tissue penetration, and cell staining for flow cytometry.

Dissociation constant Kd and affinity.