

For Research Use Only

SUMO-Tag-Trap Magnetic Agarose, Kit for Immunoprecipitation



www.ptglab.com

Catalog Number: **sutmak**

Basic Information

Catalog Number:
sutmak

Applications:
IP, Co-IP

Conjugate:
Magnetic Agarose beads; ~40 μ m (cross-linked 6% magnetic agarose beads)

Host:
Alpaca

Type:
Nanobody

Class:
Recombinant

Description

The ChromoTek SUMO-Tag-Trap Magnetic Agarose Kit consists of an anti-SUMO-tag Nanobody/VHH, which is coupled to magnetic agarose beads. It also contains lysis, wash, and elution buffers that can be used for the immunoprecipitation of SUMO-tagged proteins from cell extracts of various organisms. It can also be used in conjunction with SUMO proteases such as SenP2 for on-bead digestion of SUMO-Tag fusion proteins to release the protein of interest.

Specificity/Target

Binds specifically to all common variants of the SUMO-Tag. The SUMO-Tag is based on Small Ubiquitin-like Modifier (SUMO) proteins of a size of ca. 11 kDa, which are covalently attached to target proteins as a post-translational modification. Fusion of the SUMO-Tag to a protein of interest (POI) may increase expression and solubility of the POI. Also, the SUMO-Tag can be specifically removed by SUMO proteases such as SenP2 without leaving non-native residues behind. At least three SUMO variants are commonly used as SUMO-Tag and are all recognized by the ChromoTek SUMO-Tag-Trap: the yeast SUMO homolog SMT3, the human SUMO3 and SUMOStar, a version of SMT3 resistant to SUMO proteases. Please note that the ChromoTek SUMO-Tag-Trap will also bind non-discriminatorily to endogenous SUMO homologs such as SUMO1, SUMO2 and SUMO3 present in human cells.

Elution buffer

2x SDS-sample buffer (Lamml)

Affinity (K_D)

50 nM

Storage

Storage:

Shipped at ambient temperature. Upon receipt store at +4°C. Stable for one year. DO not freeze!

Storage Buffer:


20% ethanol

For technical support and original validation data for this product please contact


T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or E: proteintech@ptglab.com
1(312) 455-8498 (outside USA) W: www.ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data

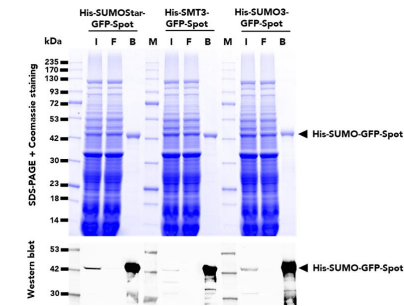

Immunoprecipitation of three different variants of SUMO-tagged, GFP fusion proteins from E. coli cell lysates using SUMO-Tag-Trap Magnetic Agarose. WB analysis was also done on samples from the Input (I), Flow-Through (F), and Bound (B) fractions using ChromoTek Mouse anti-Spot-tag Antibody (28a5) and Multi-rAb CoraLite Plus 647-Goat Anti-Mouse Recombinant Secondary Antibody (RGAM005).

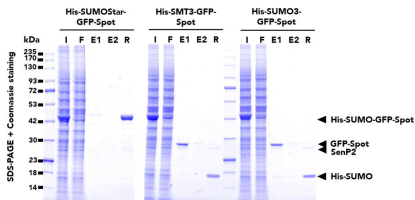

Immunoprecipitation of three different variants of SUMO-tagged, GFP fusion proteins from HEK293T cell lysates using SUMO-Tag-Trap Magnetic Agarose (sutma) followed by on-bead digest and elution with SenP2 protease. The majority of the cleaved protein is released after the 1st elution (E1), with the SenP2 and SUMO-tag remaining with beads in the residual (R) fraction. SUMOStar is resistant to SenP2 and was used as a control.